

Multi-Patient Unit



FEATURE	BENEFIT
All Metal Construction	Rugged and durable for use in multiple patient situations
Multiple Oxygen Outlets – Supports Up to Seven Patients	Six adjustable oxygen therapy outlets plus one DISS self-sealing outlet for attachment of oxygen powered resuscitators and ventilators.
25 foot Supply Hose	The unit can be placed where it is needed.
Supplied with Oxygen Therapy Masks and Tubing	Ready for use in one complete system
Soft Sided Carrying Case with Shoulder Strap	Easy to “grab and go.”

Specifications

Input pressure	50 PSI Operating	Temperature	0 to 122°F (18 to 50 °C)
Storage Temperature	-40 to 140°F (-40 to 60°C)	Storage and Operating Use	15 to 95% Humidity
Input Connection	9/16 DISS	Output Connection	9/16” DISS self-sealing outlet
Therapy Flow Connections	Therapy Barb	Therapy Flow Range	0 to 15 L/min
Weight (multi patient unit only)	19.2 oz. / 0.45 kg (approx..)	Dimensions in Inches	8.0 X 5.0 X 2.0 in. (approx.) 203 X 127 X 50 mm (approx.)




O-Two Medical Technologies Inc.
45A, Armthorpe Road, Brampton, Ontario, Canada, L6T 5M4
Office: +1 905 792 OTWO (6896)
N.A. Toll Free: +1 800 387 3405 Fax: +1 905 799 1339
Email: resuscitation@otwo.com Web: otwo.com

Multi-Patient Unit

Designed for 6 patients needing oxygen therapy simultaneously:

9/16" DISS inlet

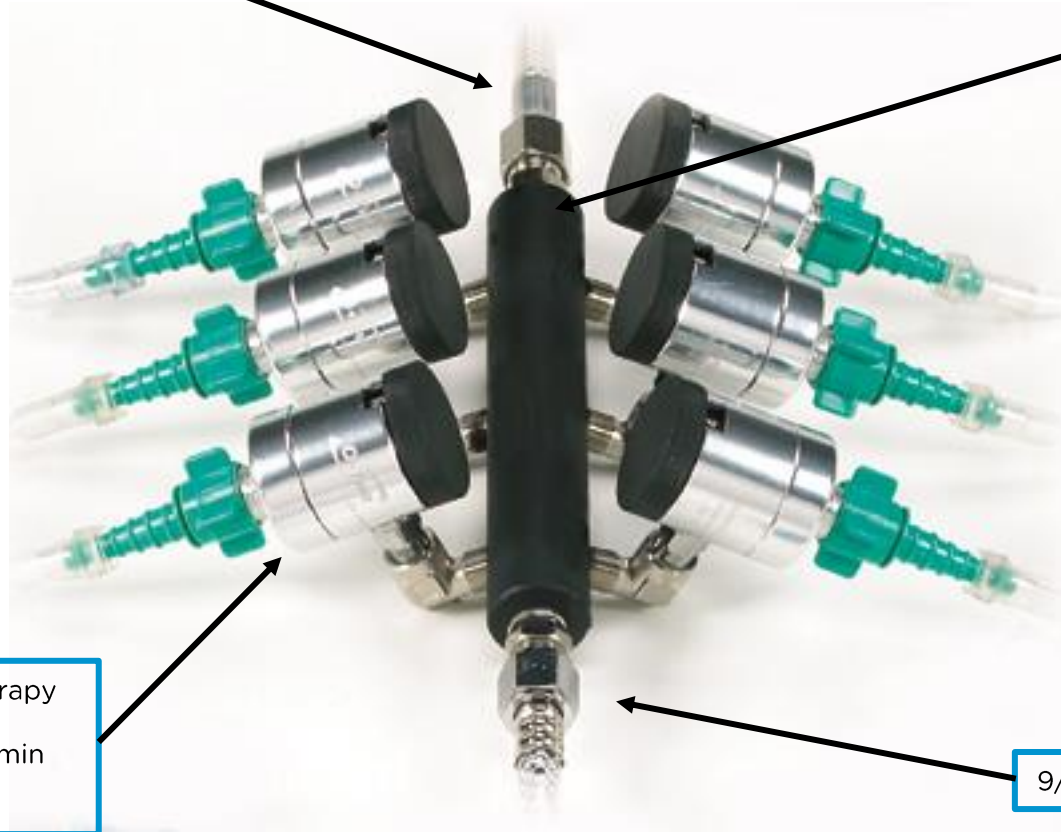

All Metal Construction

- Provides Uniform Pressure to All Flow Controllers

Six individually adjustable therapy flow controls

- Each provide 0 - 15 L/min flow rate

9/16" DISS self-sealing outlet

Kit Includes:

- A twenty-five foot supply hose
- Six adult medium concentration therapy facemasks with sure-flow tubing
- Soft-pack carrying case with shoulder strap